


Water-heater



www.horse-dental-equipment.com

Water-heater

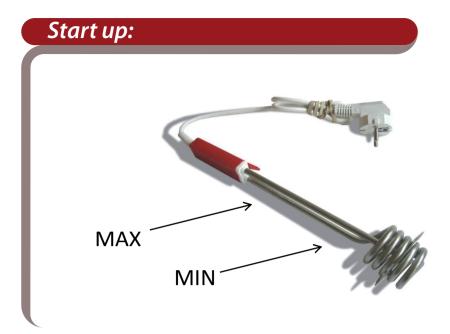

Only use metal or glass containers with a minimum capacity of 100ml (with a minimum diameter of 70mm).

Fill your container with the liquid you wish to heat.

Partially immerse the water heater so that the liquid level is between the MIN and MAX marks on the heater

Plug in the socket

As soon as the water boils or arrives at the desired temperature, first disconnect the plug and then remove the heater from the water.

After being unplugged the heater stays hot. Do not touch it or put it on fragile or combustible surfaces.

Recommandations:

This device should not be left unattended when in use.

Do not use an external timer or a remote control system.

Risks if you leave it unattended: the liquid will evaporate due to the heat. After a few minutes of overheating, the heater will be put permanently out of use by the thermal safety in the device. In order to continue use your heater for a long time, make sure that the liquid level never falls below the MINI mark.

Maintenance:



Disconnect and clean the appliance after use with a soft cloth.

